



FLOCKTAGGER APPLICATOR




LOADING THE APPLICATOR

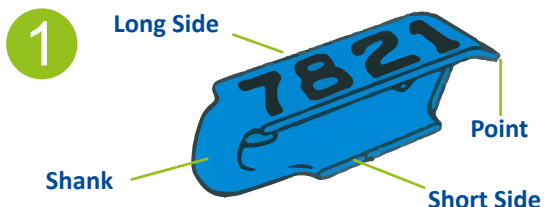

1 Load the tag with the point towards the applicator. The short side locates in the loading slot.

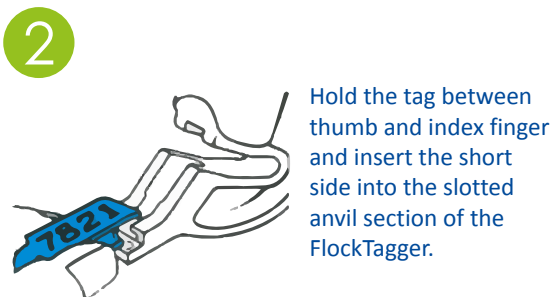

APPLYING THE TAG

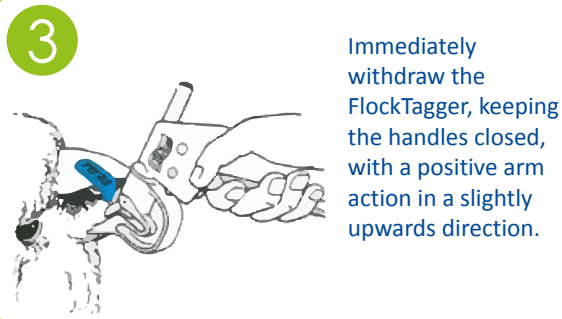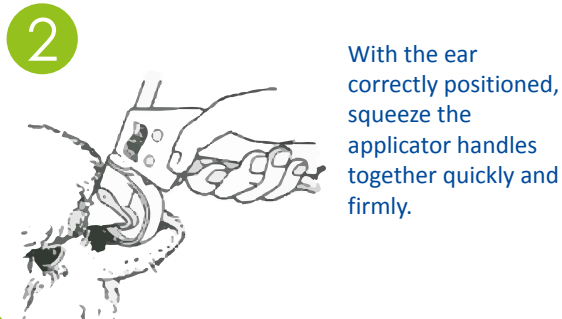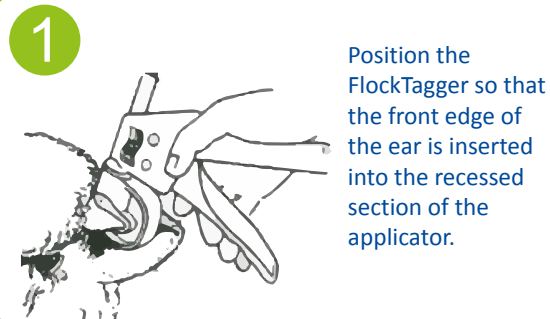

IMPORTANT:

Squeeze the handles together quickly and firmly, keeping the handles together withdraw the FlockTagger with a positive and prompt movement of the arm. Avoid pulling sideways. It is important to always withdraw the applicator in the correct direction, a little practice may be required to provide the correct technique.

Refer to the FlockTag placement diagram for correct tag positioning on the ear.



TAG POSITIONION:

Insert the tag in from the front of the top edge of the ear, to the full depth of the applicator or less depending upon the age or breed of animal. Always allow for future growth and or differing ear sizes.

Make sure the tag is far enough in from the ear tip to prevent it from swinging around past the tip.



Contact our friendly Customer Service Team

1300 138 247

www.allflex.com.au



Identify the Difference